



Cobalt Oxide Nanoparticles COON 1601

Description:

Co_3O_4 has a cubic spinel crystal structure in which the Co^{2+} ions occupy the tetrahedral sites and the Co^{3+} ions the octahedral sites. The Co^{3+} ions at the octahedral sites are diamagnetic in the octahedral crystal field. The Co^{2+} ions at the tetrahedral sites form an antiferromagnetic sublattice with a diamond structure. Co_3O_4 spinel, have been intensively studied in photocatalytic systems because of their superior properties, good stability, and high abundance. Co_3O_4 nanoparticles exhibit weak ferromagnetic behavior and CoO nanocrystals display super paramagnetism or weak ferromagnetism, whereas bulk CoO is antiferromagnetic. Co_3O_4 is a magnetic p-type semiconductor.

Characterization	
CAS	1308-06-7
Stock No.	COON 1601
Molecular formula	Co_3O_4
Molecular weight (g/mol)	240.80
Form	Powder
Color	Black
Morphology	
Crystal structure	FCC
Size range (nm)	
Melting Point	895°C
Boiling Point	900°C
Density (g/cm ³)	6.11
Solubility	Insoluble



Image of Cobalt oxide nanopowder (COON 1601)


Note: product specifications are subject to amendment and may change over time.

Applications (but not limited to the following):


The cobalt spinel compounds can act as efficient catalysts in a lot of heterogeneous chemical processes. Nanoparticles of Co_3O_4 are promising materials for electronic devices, gas sensors, magnetic materials, electrochromic devices, electrochemical systems, and high-temperature solar selective absorbers. CoO also shows interesting properties and has applications as gas sensors and as anodes of lithium-ion batteries.

Safety:

Avoid breathing dust.
Always use protective gloves and safety glasses.
Wash with soap and water after exposure.
Refer to MSDS prior to handling this material.

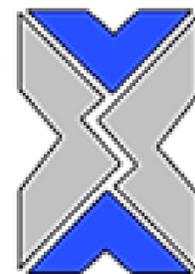

Cobalt Oxide Nanoparticles COON 1601

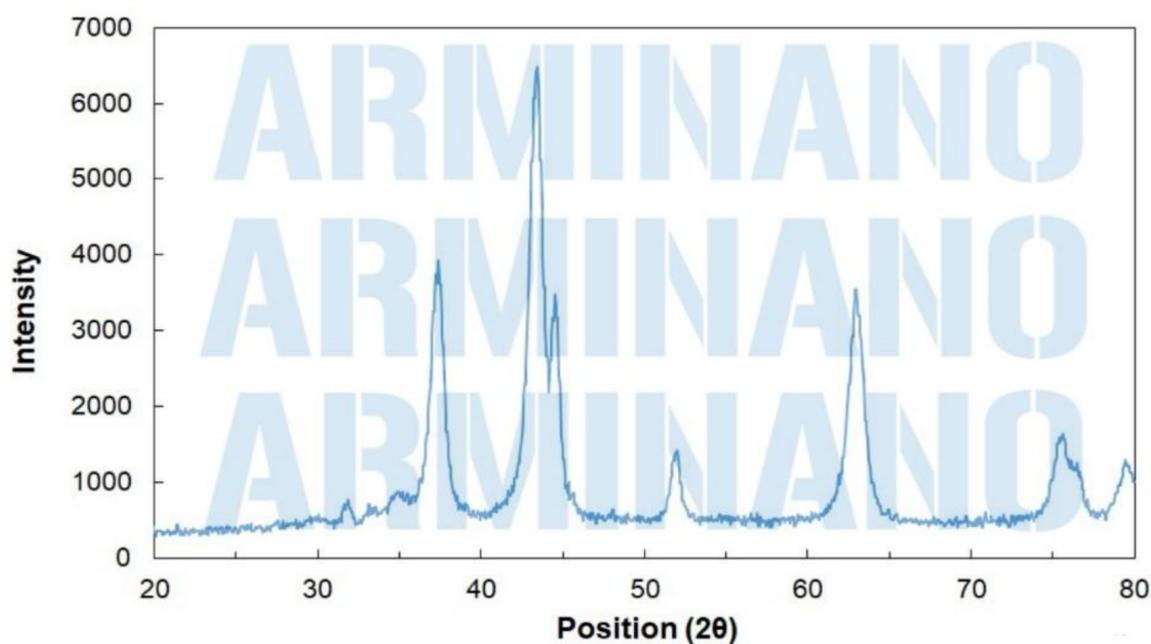

XRD pattern of NOP1201

Storage:

- Keep it in cool dry place and tightly closed container.
- Ensure good ventilation at the workplace.
- The product is not flammable.
- Avoid direct sunlight.
- Store away oxidizing agents.
- To disperse nanoparticles sonication could be used.

Shelf life:

When stored as specified the product is stable for at least 6 months.

